# **NF 048**

# MICROLITRE AND HAEMATOCRIT CENTRIFUGE




NF 048 is a compact centrifuge that can be used both as microlitre and haematocrit centrifuge.

With its microlitre rotor accepting small capacity tubes up to 2 ml, NF 048 is ideal for hospitals and research laboratories and for applications such as pelleting of DNA and protein, DNA amplification, enzyme tests, centrifugation of cells, yeast and microorganisms at high speed.

When NF 048 is operated with its haematocrit rotor, it can be used for determining the haematocrit volume by means of the centrifugation of blood samples in capillary tubes.

The body, lid and chamber of NF 048 is made of epoxy-polyester powder coated steel which means that the centrifuge is resistant to external and internal effects.

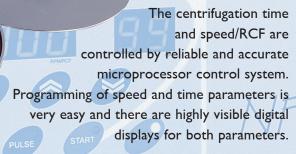

NF 048 ensures minimum increase of sample temperature by means of the continuous air flow system through the air channels in the lid.

By means of the lid locking system, the lid can not be opened while the rotor is spinning and the program can not be started if the lid is open or not closed well.

Quiet and brushless induction motor of NF 048 has a very high performance and the motor accelerates the rotor to the maximum speed in a very short time. The motor is protected against over heating. There is a pulse key for short centrifugation requirements.

The microlitre rotor is made of polypropylene which is resistant to a wide range of chemicals and can be autoclaved.





# **TECHNICAL SPECIFICATIONS**

|                               | NF 048  |  |  |
|-------------------------------|---|--|--|
| Maximum Speed                 | Angle (Microlitre): 14,000 rpm Haematocrit : 12,000 rpm         |  |  |
| Maximum RCF                   | Angle (Microlitre): 18,188xg Haematocrit : 14,811xg             |  |  |
| Maximum Capacity              | Angle (Microlitre): 24x1.5/2 ml Haematocrit : 24xcapillary tube |  |  |
| Control System                | Programmable Microprocessor Control System                      |  |  |
| Speed Set Range               | 1,000 - 14,000 rpm  |  |  |
| Speed Set Step                | 10 rpm  |  |  |
| Timer Set Range               | I - 99 min. and hold position                                   |  |  |
| Timer Set Step                | I minute  |  |  |
| Motor                         | Maintenance-free, brushless induction motor                     |  |  |
| Power Supply                  | 230 V, 50/60 Hz.  |  |  |
| Maximum Power                 | 450 W   |  |  |
| Overall Dimensions (WxDxH) mm | 275x355x240   |  |  |
| Packing Dimensions (WxDxH) mm | 350x420x460   |  |  |
| Net/Gross Weight, kg          | 13/17   |  |  |
|                               |   |  |  |

#### **ROTORS**

| Cat. No  | Rotor Type        | Capacity          | Max. Tube Dia. (mm) | Radius (mm) | Max.Speed (rpm) | Max. RCFxg |
|----------|-------------------|-------------------|---------------------|-------------|-----------------|------------|
| B 50 019 | Angle(Microlitre) | 24x1.5/2 ml       | П                   | 83          | 14,000          | 18,188     |
| B 50 020 | Haematocrit       | 24xcapillary tube | 0.75                | 92          | 12,000          | 14,811     |

### **ADAPTORS FOR 24x1.5/2 ml ANGLE (MICROLITRE) ROTOR**

| Cat. No. | Description | Capacity         | Max. Tube Dia. (mm) |
|----------|-------------|------------------|---------------------|
| A 14 005 | Adaptor     | l×500/800 μl     | 8                   |
| A 14 006 | Adaptor     | Ix200 μIPCR      | 6.5                 |
| A 14 007 | Adaptor     | I×250/400/700 μI | 6                   |





NF 048



## NÜVE SANAYİ MALZEMELERİ İMALAT VE TİCARET A.Ş.

Esenboğa Yolu, 22 km. Akyurt 06750 ANKARA TURKEY Tel : (90.312) 399 28 30 Fax : (90.312) 399 21 97 www.nuve.com.tr sales@nuve.com.tr

ISO 9001:2008 CE



We reserve the right to modify the design and description